

On scattering effects for volume sources in low-energy photon spectrometry

Marie-Christine Lépy, Laurine Brondeau, Laurent Ferreux, Sylvie Pierre

▶ To cite this version:

Marie-Christine Lépy, Laurine Brondeau, Laurent Ferreux, Sylvie Pierre. On scattering effects for volume sources in low-energy photon spectrometry. Applied Radiation and Isotopes, 2013, 81, pp.71 - 75. 10.1016/j.apradiso.2013.03.071. cea-01802412

HAL Id: cea-01802412 https://cea.hal.science/cea-01802412

Submitted on 24 Mar 2020

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

ON SCATTERING EFFECTS FOR VOLUME SOURCES IN LOW-ENERGY PHOTON SPECTROMETRY

Marie-Christine Lépy, Laurine Brondeau, Laurent Ferreux, Sylvie Pierre CEA, LIST, Laboratoire National Henri Becquerel, Gif-sur-Yvette, F-91191, France

ABSTRACT

In this study, different aspects of the Compton scattering inside volume sources are illustrated using experimental approach and Monte Carlo simulation. For the low-energy range (below 100 keV) scattered events represents around 30 % of the whole spectrum. Influence of the source-detector geometry is discussed. The scattering effects induce strong differences in spectrum shape for different geometrical conditions. This should influence efficiency transfer factors. A new approach is proposed, including the scattered events, to avoid complex peak area determination.

Keywords: Compton scattering, low-energy, Monte Carlo simulation, PENELOPE, efficiency calibration

Highlights:

Compton scattering is unavoidable effect in volume sources.

In the low-energy range, this induces registration of scattered photons close to the full-energy peak.

The respective contributions of the source material and container are studied for different geometrical conditions.

The consequence on the full-energy peak area determination is discussed.

INTRODUCTION

Gamma-ray and X-ray spectrometry in the low-energy range (20-150 keV) suffers a number of difficulties that increase uncertainties associated with quantitative results derived from this region. The main point is the full-energy peak (FEP) efficiency calibration in this energy range, where the curvature of the efficiency versus energy function changes, around its maximum value. Nowadays, Monte Carlo simulation is currently used to calculate efficiencies; however, the validity of the calculation must be carefully checked by comparing the calculated efficiencies with experimental values. Moreover, photon spectra are strongly distorted by scattering effects. Compton scattering occurs either in the source (matrix and container for a volume source), in the source holder, in the detector itself or its shielding. The resulting scattered « peak » is a combination of the relevant scattering angles, thus giving a continuous wide bump close to the FEP: the lower the energy, the closer the scattered events. Consequently, this effect prevents accurate determination of the FEP area that is used either for the calibration or for quantitative measurements. In this study, qualitative information of the scattering effect is derived from Monte Carlo simulation and illustrated with experimental approach.

1. Experimental

This study concentrates on spectra recorded with a coaxial n-type high purity germanium (HPGe) detector, linked to analogue electronics, described in (Ferreux *et al.*, 2012). The detector has an active volume of 100 cm³ and a beryllium window with a thickness of 500 µm. Three types of sources with an activity of about 40 kBq are used:

• Point sources, constituted by a weighted drop of standard radioactive solution included between one Mylar® foil adhesive tape (top), these foils being maintained by a plastic ring;

- Liquid volume sources in plastic cylindrical containers, with volumes 50 cm³ and 500 cm³;
- Gaseous volume sources in a stainless steel 100 cm³ cylindrical container.

1.1 Efficiency calibration in the low-energy range

To get experimental efficiency calibration, one of the difficulties is the lack of reliable photon emission references, since standard radionuclides emitting photons at energies below 100 keV are scarce. The Laboratoire National Henri Becquerel (LNHB) develops and maintains the "Table des Radionucléides" (Bé et *al.*, 2011) providing carefully evaluated radionuclide decay data with extensive information on the data origin. For the energy range of interest, these include a few gamma emitters (57Co, 75Se, 109Cd, 133Ba, 210Pb, 241Am) with some questionable photons emission intensities. In addition, KX-rays of these radionuclides and some others can be used, but it must be noted that their emission intensities are derived from computation, using K capture or conversion-electron probabilities and fluorescence yields. For this last parameter, strong correlation exists in the database, inducing correlated X-ray emission intensities values.

The efficiency calibration of the detector is obtained using standard solutions and is accurately established for point and 50 cm³ liquid sources, at 10 cm from the detector window, and for 500 cm³ liquid sources at 8.33 cm.

1.2 Scattering effect in the low-energy range

A previous study (Plagnard *et al.*, 2008) was carried out to determine the contribution of surroundings of point sources to the scattered events. The results indicate that precautions

should be taken to reduce unwanted bumps, particularly by reducing the material immediately around the radioactive deposit. Furthermore the use of a source collimator significantly reduces the scattering effect in the resulting spectra. However, these favourable geometrical conditions cannot always be attained, depending on the sample shape, source-to-detector distance, etc. Particularly, in low-level measurements, volume samples are currently used, thus scattering in the sample volume is unavoidable and must be considered. Figure 1 compares spectra obtained with ²⁴¹Am in the 5 to 70 keV energy range, using a point source (black) and a 50 cm³ volume source (red) at 10 cm from the detector window. On the left of the 59.6 keV peak, clearly, for the point source, a scattering bump appears around 52.5 keV, whereas for the volume source, the scattering effect gives a continuous distribution ranging from 40 keV up to the full-energy peak (FEP) energy.

Spectra obtained with ¹³³Ba in the same conditions are presented in Figure 2, in the 10 to 90 keV energy range where are both the caesium K X-rays and the 79-80 keV doublet of ¹³³Ba. Here the features seen in ²⁴¹Am spectra are also shown for both energy regions. The scattered events induce a relatively distinct bump for the point source, centred at 28.8 and around 68 keV for the 30 and 80 keV, respectively. For the volume source, they lead to a continuous distribution, starting from around 27 and 60 keV up-to the FEP energy, thus strongly disturbing the FEP background.

2. Monte Carlo simulations

To estimate the respective contributions of the source and container materials in the scattering distributions, Monte Carlo simulations were performed, using the PENELOPE code (Salvat *et al.*, 2009). PENELOPE2008 performs simulation of coupled electron and photon transport in arbitrary materials; it includes several generic main programs allowing an easy

implementation. The PENMAIN routine performs particle tracking in material systems consisting of homogeneous regions (bodies) limited by quadric surfaces. One of the output files provides the distribution of absorbed energy in selected bodies (energy deposition spectrum). The code is written in FORTRAN and has been adapted to track any Compton scattering event in elements surrounding the source before energy deposition in the germanium crystal. The main sites of scattered events are then identified and the resulting energy distributions are given for each material.

2.1 Volume and distance effects

Figures 3-5 present the results of the Monte Carlo simulation for 60 keV photons, for three geometries under study: point source, 50 and 500 cm³ volume (water) sources. In the last ones, the events consecutive to scattering, either in the liquid or in the container, are plotted. These simulated spectra can be compared with the americium-241 spectra (Figure 1). For 60 keV photons, the relative contribution of Compton scattering (total spectrum minus full-energy peak) represents 11.3 % of whole spectrum for the point source, whereas it amounts to 32.5 % and 45.3 % for 50 and 500 cm³ volume sources, respectively. It can also be stated that, for the point source, the scattered events lead to a peaked distribution, centred around 52 keV, far from the FEP. For the volume sources, this distribution starts around 40 keV up to the incident photons energy. In addition, the detailed analysis demonstrates that the scattering effect in the liquid is the main contribution to the continuous bump obtained in experimental spectra, whereas scattering in the container has a less significant impact.

Figure 6 presents the simulation results for 60 keV photons emitted from a 50 cm³ volume source at the contact of window, what can be compared with Figure 4 where the source is at 10 cm. Due to the shorter source-to-detector distance, the ratio of scattered events in the

whole spectrum represents about 28 % when the container is at the contact, and amounts up to 32.5 % when it is at 10 cm from the detector window.

2.2 Energy effect

For the 50 cm³ volume source at 10 cm from the detector window, the Monte Carlo simulation has been performed with different incident energies. It was stated that the relative contribution of scattering increases with energy, representing 30.3, 32.5, 35.1 and 37.2 % of whole spectrum, for energies 40, 60, 80 and 100 keV, respectively.

2.1 Material effect

In the following, the study concentrates on scattering in the source or its container, which corresponds to the high energy part of the spectra, on the left side of the FEP peak. The influence of the source material is examined for 60 keV photons either coming from water, silica or nitrogen included in a 50 cm³ plastic container at 10 cm from the detector window. The respective results are presented in Figures 4,7 and 8. For water and silica, the scattering distribution is mainly due to interactions in the matrix source, whereas for the gas source, scattering in the container is the dominant effect.

Finally, the influence of the container material has been examined by comparing spectra obtained with 80 keV photons emitted by water either in plastic or aluminium 50 cm³ container. The relative contribution of Compton scattering represents about 34 % and 37 % of whole spectrum for the plastic and aluminium containers, respectively. Scattering in the water is of the same size of order (about 22.5 %) for both cases. However, it is twice higher in the aluminium container (about 6 %) than in the plastic one.

3. Application to quantitative results

3.1 Influence on FEP area determination

In the low-energy region, scattering effect strongly disturbs spectrum shape and the resulting high background prevents from accurately determining the FEP areas. In the traditional approach to determine FEP, background subtraction is performed using either linear or step function. The spectrum presented in Figure 2 (133 Ba – Volume source) can be processed according to this approach, using COLEGRAM software (Ruellan *et al.*, 1996). However, according to the above-mentioned results, scattering effect should be taken into account and the processing should include this feature. On the quantitative aspect, the resulting peak areas (sum of 79.6 and 81.0 keV) give respectively 1.741 x 10^6 ("normal" background subtraction) and 1.786 x 10^6 (background subtraction including scattering), which gives a relative difference of more than 2.5 %.

3.2 Efficiency transfer

Monte Carlo simulation can be used to compute efficiency for different geometrical conditions, thus allowing derivation of transfer factors that can be used to determine efficiency values for a measurement geometry different from the calibration one (Vidmar *et al.*, 2010). However, as discussed above, it can be difficult to achieve accurate peak area determination due to the scattering contribution that strongly depends on the volume geometry. Such a difference between calibration and measurement geometries is found in the case of ¹³³Ba and ¹³³Xe: both decays towards excited levels of ¹³³Cs thus emitting the same photon energies. Consequently, efficiency calibration obtained from ¹³³Ba should be easily transferred to measure ¹³³Xe activity. As a practical application, it was required to determine

the activity of a ¹³³Xe sample included in a stainless steel container placed at 6.8 cm at from the detector window. Since no calibration is available for this geometry, the activity measurement must include an efficiency transfer factor from ¹³³Ba calibration obtained with a 50 cm³ liquid sources. But, due to the nature of the elements (barium is included in a solution, whereas xenon is a gas), the geometrical conditions are different and so are relevant to the experimental spectra, with very strong scattering in the case of xenon.

Another approach can thus be suggested for such volume samples: the transfer factor should consider not only interactions inducing deposit of the whole initial energy of the photons (FEP), but also those following scattering out of the detector. This is illustrated in Figure 9-a, that presents Monte Carlo simulation results for 81 keV photons following the decay of ¹³³Xe (gas in a 100 cm³ stainless steel container). The energy region between 37 and 81 keV allow calculation of a "transfer factor including scattering" for the 81 keV peak. This can be applied to the experimental spectrum (Figure 9-b) using the net area of the energy region plotted in red. However, this approach is valid in such a specific case where the full-energy peak is well isolated with no other contribution in the scattering region.

4. Conclusion

In this study, different aspects of the Compton scattering inside volume sources are illustrated using experimental approach and Monte Carlo simulation. For the low-energy range (below 100 keV), the scattering contribution is significant in gamma-ray spectra and represents around 30 % of the recorded events. Influence of the source-detector geometry is discussed: for liquid sources, scattering in the source are dominant compared to scattering in the container. On the contrary for gas sources, scattering in the container is the major

contribution. These scattering effects induce strong differences in spectrum shape for different geometrical conditions. This should influence efficiency transfer factors and a new approach is proposed which is valid in specific conditions where the full-energy peak and associated scattered distribution are well isolated. In these cases, the efficiency transfer factor can be computed including the scattered events, what avoids complex peak area determination to get accurate quantitative results.

Figures:

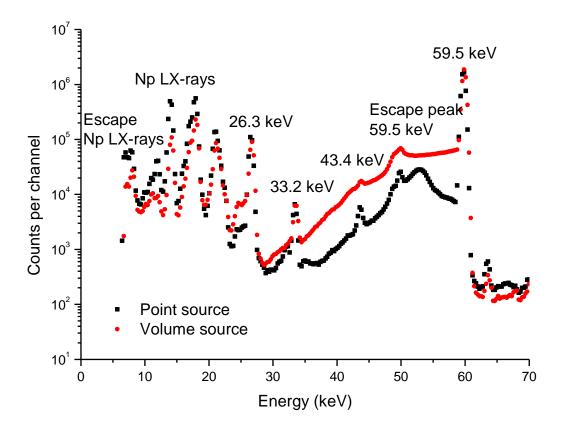


Figure 1: Spectra recorded with ²⁴¹Am at 10 cm from the detector window: Black: Point source; Red: Water in a 50 cm³ plastic container

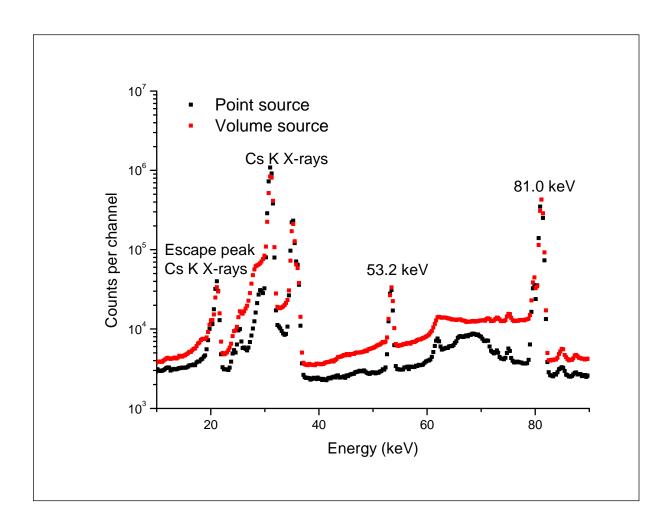


Figure 2: Spectra recorded with ¹³³Ba at 10 cm from the detector window: Black: Point source – Red: Water in a 50 cm³ plastic container

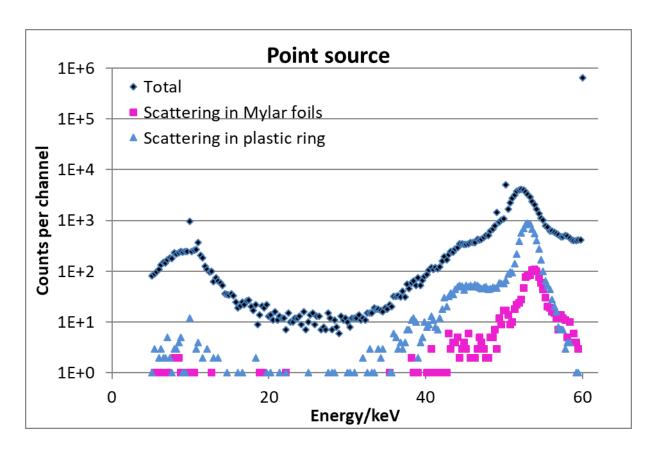


Figure 3: Monte Carlo simulation for 60 keV photons emitted by a point source at 10 cm from the detector window

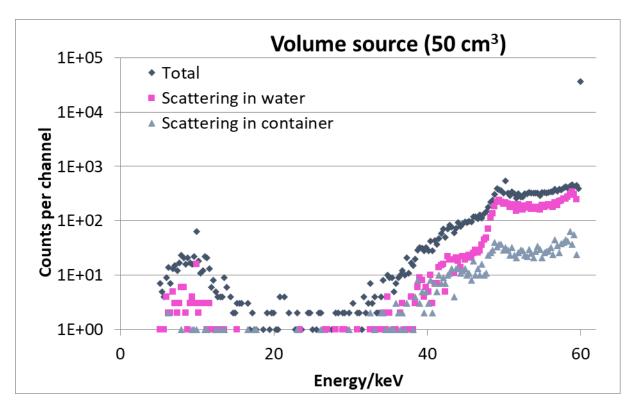


Figure 4: Monte Carlo simulation for 60 keV photons emitted by water in a 50 cm³ plastic container at 10 cm from the detector window

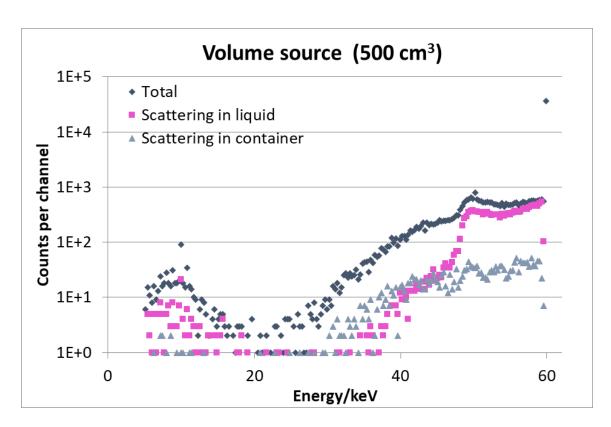


Figure 5: Monte Carlo simulation for 60 keV photons emitted by water in a 500 cm³ plastic container at 10 cm from the detector window

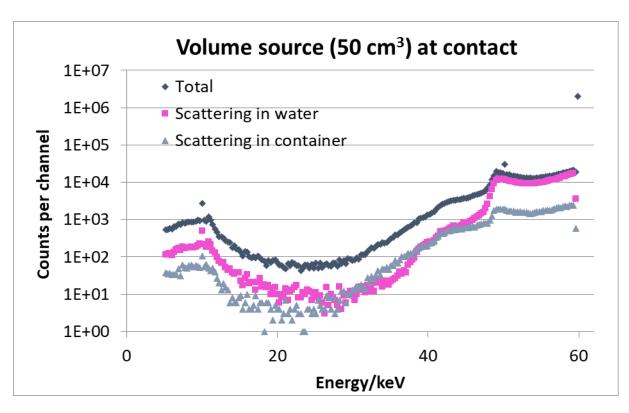


Figure 6: Monte Carlo simulation for 60 keV photons emitted by water in a 50 cm³ plastic container at contact of the detector window

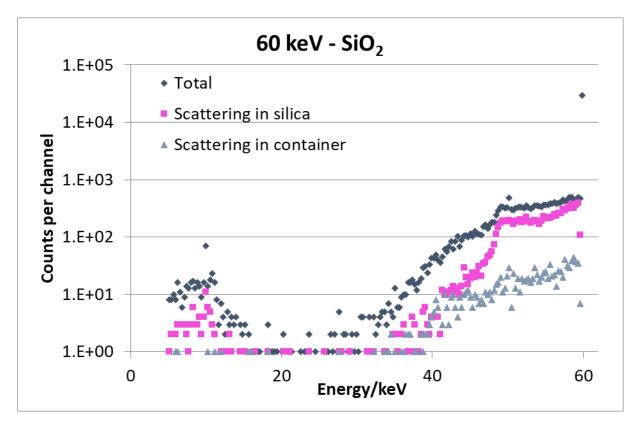


Figure 7: Monte Carlo simulation for 60 keV photons emitted by silica in a 50 cm³ plastic container, at 10 cm from the detector window

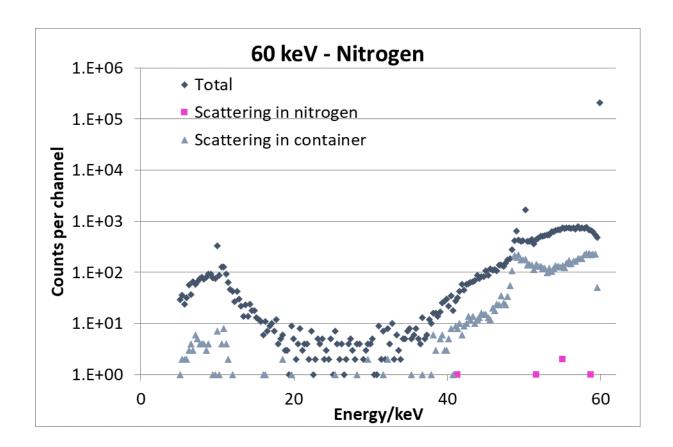


Figure 8: Monte Carlo simulation for 60 keV photons emitted by nitrogen in a 50 cm³ plastic container, at 10 cm from the detector window

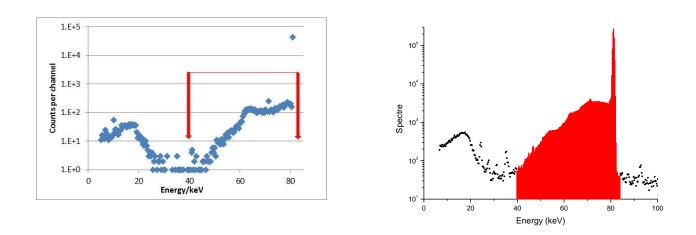


Figure 9: ¹³³Xe: - a: Monte Carlo simulation of 80 keV photons interaction – b: Experimental spectrum

References

Bé, M.-M., Chisté, V., Dulieu, Mougeot, X., C., Browne, E., Chechev, V., Kuzmenko, N., Kondev, F., Nichols, A., Luca, A., Galan, M., Arinc, A., Huang, X., 2011. Table of Radionuclides, Monographie BIPM-5, ISBN 13 978 92 822 2234 8 (set) and ISBN 13 978 92 822 2235 5 (CD), CEA/LNE-LNHB, 91191 Gif-sur-Yvette, France and BIPM, Pavillon de Breteuil, 92312 Sèvres, France.

Ferreux, L., Pierre, S., Thanh, T.T, Lépy, M.-C., 2012. Validation of efficiency transfer for Marinelly geometries. These proceedings

Plagnard, J. Hamon, C., Lépy, M.-C., 2008. Study of scattering effects in gamma-ray spectrometry. Applied Radiation and Isotopes 66, 769-773.

Ruellan, H., Lépy, M.-C., Etcheverry, M., Morel, J., 1996. A new spectra processing code applied to the analysis of ²³⁵U and ²³⁸U in the 60 to 200 keV energy range. Nucl. Instrum. Meth. in Phys. Res. A369, 651-656.

Salvat, F., Fernandez-Varea, J.M., Sempau, J., 2009. PENELOPE-2008: A code system for Monte Carlo Simulation of Electron and Photon transport, OECD Nuclear Energy Agency, Issy-les-moulineaux, France.

Vidmar, T., Çelik, N., Cornejo Díaz, N., Dlabac, A., Ewa, I.O.B., Carrazana González, J.A., Hult, M., Jovanović, S., Lépy, M-C., Mihaljević, N., Sima, O., Tzika, F., Jurado Vargas, M.,

Vasilopoulou, T. and Vidmar, G. 2010. Testing efficiency transfer codes for equivalence. Applied Radiation and Isotopes 68, 355-359.